



ALVERNO LABORATORIES

SPECIMEN COLLECTION GUIDE MOLECULAR

Specimen collection instructions and guidelines may change. We recommend checking the collection manual for the most up-to-date collection information at www.AlvernoLabs.com

Reviewed December 2020

Reach us at 800-937-5521

THIN PREP PAP



ORDER# 144598

⇒ HPV

⇒ GC/Chlamydia

ABBOTT MULTI-COLLECT SPECIMEN COLLECTION KIT



ORDER# 176526

Intended for the collection and transportation of male and female, swab and urine specimens for the detection of Chlamydia trachomatis and Neisseria gonorrhoeae

BBL CULTURE SWAB LIQUID STUART

ORDER# 193216

OR E-SWAB

ORDER# 132433

⇒ MRSA or MSSA by PCR



BD AFFIRM VAGINITIS PANEL

ORDER# 108401

⇒ Gardnerella

⇒ Trichomonas

⇒ Candida sp.*

*does NOT identify the species of Candida yeast




BD Max Vaginitis Panel

ORDER# 181279

⇒ Bacterial Vaginitis

⇒ Trichomonas

⇒ Candida species will be identified

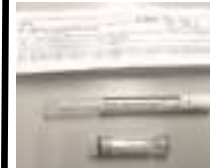

E-SWAB

ORDER# 132433

⇒ Group B Strep by PCR

⇒ Vaginal/Rectal E-SWAB



STERILE SCREW CAP CUP (no additive)

ORDER# 153706

⇒ Respiratory Wash/Aspirate

* Flu A / Flu B / RSV panel

* Flu A / Flu B panel

* Respiratory 21 Pathogen Panel

* B. pertussis and parapertussis

* Mycoplasma pneumoniae

⇒ HSV 1 and 2

* Respiratory Wash / CSF

⇒ Varicella Zoster Virus (VZV)

* CSF

⇒ C. diff by PCR

* Stool

⇒ Enteric Pathogen Panels

* Stool (bacterial and parasitic)

⇒ Norovirus

* Stool

⇒ CMV Qualitative by PCR:

* CSF, Urine, Respiratory Wash



BD UTM (Universal Transport Media)

ORDER# 111005

By PCR:

⇒ Varicella zoster virus (VZV)

⇒ Herpes virus 1&2

⇒ Flu A / Flu B / RSV virus panel

⇒ Flu A / Flu B virus panel

⇒ Mycoplasma pneumoniae

⇒ B pertussis and parapertussis

⇒ Respiratory 21 Pathogen Panel

⇒ SARS-CoV-2



LAVENDER TOP TUBE: WHOLE BLOOD

ORDER# 193813

⇒ Factor V Leiden Mutation by Real Time PCR

⇒ Factor II Mutation by Real Time PCR

⇒ Cystic Fibrosis screen: original tube only

⇒ CD4/CD8 (T-Cell Subsets)

⇒ Lymphocyte Subsets (T, B, and NK cells)

⇒ Lymphoma/Leukemia Immunophenotyping by FlowCytometry

⇒ Fetomaternal Bleed by FlowCytometry

⇒ HLA B27 Detection by FlowCytometry



UROVYSION

ORDER# 135884

⇒ For bladder cancer test collection use Urocyte Kit



PARA-PAK SAF FIXATIVE

ORDER# 11939

⇒ Enteric Parasite Panel

PARA-PAK C&S CARY-BLAIR MEDIA

ORDER# 23205

⇒ Gastrointestinal 22 Pathogen Panel

⇒ Enteric Bacterial Panel

⇒ Extended Enteric Bacterial Panel

⇒ Norovirus



LAVENDER TOP TUBE: PLASMA*

ORDER# 193813

⇒ CMV Qualitative PCR

⇒ CMV Quantitative PCR

⇒ HCV Quantitative by PCR

⇒ HIV Quantitative by PCR**

⇒ HBV Quantitative By PCR

** requires plasma from 2 lav tubes

*Submit plasma in a transfer tube



GREEN TOP TUBE: SODIUM HEPARIN

ORDER# 193842

⇒ Hematological FISH testing for peripheral blood or bone marrow aspirate

